



User Guide for GCE Zen-O™ and Zen-O™ Lite External Battery Charger

Overview

This guide provides instructions on how to use the GCE Zen-O™ and Zen-O™ Lite External Battery Charger effectively, based on the information provided on the unit itself and standard battery charging practices.

Contents

1. Charging the Battery
2. LED Indicator Status
3. Important Safety Notes
4. Troubleshooting

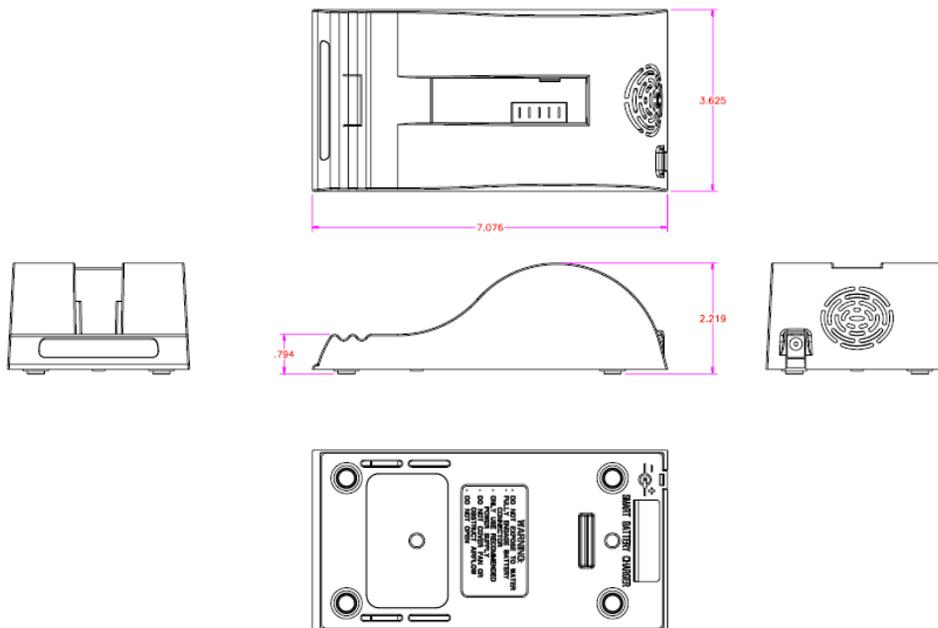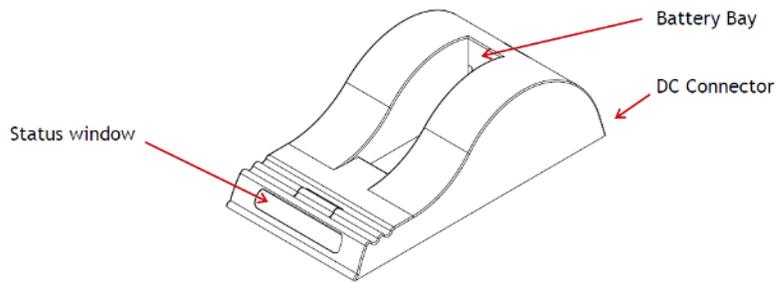

1. Charging the Battery

Step-by-Step Instructions:

1. **Prepare the Charger:**
 - Ensure the charger is placed on a flat, stable surface.
 - Connect the charger to the recommended power supply (AC 110-240V).
2. **Insert the Battery:**
 - Align the battery with the charger's connector.
 - Push the battery gently into the charger until it clicks into place, ensuring a secure connection.
3. **Begin Charging:**
 - Once the battery is connected, the charger will automatically begin the charging process.

2. LED Indicator Status

The charger features an LED indicator that provides information about the charging status:

- **Green Flashing Light:** The charger is in the process of charging the battery.
- **Green Solid Light:** The battery is fully charged and ready for use.
- **Red Solid Light:** There is a fault with the battery.

3. Important Safety Notes

- **Do not expose** the charger or power supply to water or liquids.
- **Do not open** the charger or power supply case; there are no user-serviceable parts inside.
- **Do not cover** the air vent or obstruct the airflow to prevent overheating.
- **Use only** the manufacturer's power supply and observe terminal polarity.
- **Place** the charger in a cool spot, away from external heat sources.

4. Troubleshooting

- If the **LED indicator shows a red solid light**, this indicates a battery fault. Try re-seating the battery or using a different battery.
- Ensure the **charger and battery connections** are clean and free of debris.
- Make sure the **power supply** is properly connected and functioning.

Using the Plastic Spacers

Plastic spacers are supplied to assist in the correct insertion of all different sizes of batteries compatible. To insert a spacer, place the bottom of the spacer in the rear corner of the battery bay. Place a thumb on the top of the spacer and push it firmly towards the back of the charger. Correct insertion should be accompanied by a click as the guides lock into place.